
# MediPendant™

# 10 STEP QUICK START GUIDE

Model #MED01


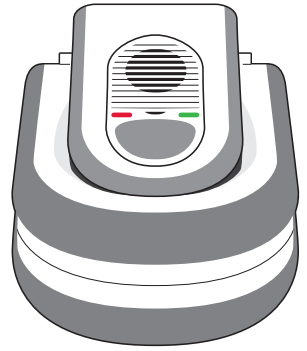
## INSTALLING YOUR MediPendant™

This is for basic installation of your **MediPendant™**. *Please see OWNER'S MANUAL for full instructions, product details, warranties and disclaimers.*

### STEP 1

Place the **MediPendant™** base station on a stable surface located near a power outlet and an available phone jack. For best results, locate the base station in an area away from household appliances which could present electrical interference. See Fig. 1

If you are using more than one **MediPendant™** system in your home, we recommend that you place the base stations at least 15 ft. apart from each other.
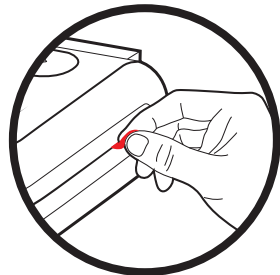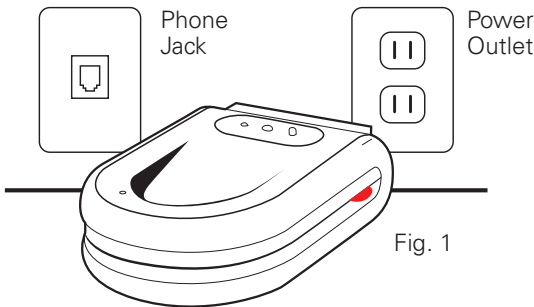
### STEP 2

Grab RED TAB located on side of base station and pull it straight out. Removing tab activates back-up battery. See Fig. 2

Unit will beep continuously and the POWER LED and PHONE LINE INDICATOR LED will both flash RED until AC/DC adapter is plugged in. See **STEP 3**.


### STEP 3

Connect the AC/DC Adapter into the **BLACK** plug in the back of the base station marked DC. See Fig. 3

### STEP 4

Connect the other end of the AC/DC Adapter into a standard power outlet. See Fig. 4

Once plugged in, your base unit will continue to beep and the RED PHONE LINE INDICATOR LED will continue to flash. The POWER LED will change from flashing RED to solid GREEN.

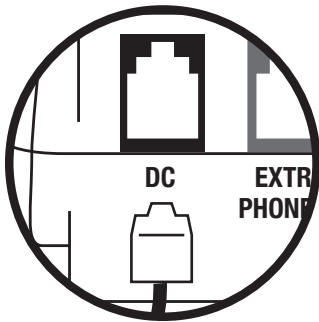STEP 3

Fig. 3

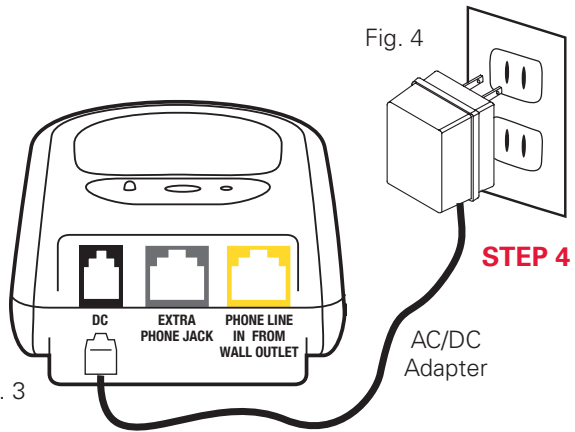
Fig. 4

STEP 4

AC/DC  
Adapter

**NOTE: DO NOT PLUG THE AC/DC ADAPTER INTO A POWER OUTLET WHICH IS CONTROLLED BY A WALL SWITCH.** If the switch is accidentally turned off, the system will begin to beep indicating a loss of power.

**NOTE:** With loss of power, the POWER LED located in the center on the top of the base station will illuminate and flash RED and the unit will beep. **TO SILENCE ANY TONE FROM THE BASE STATION,** simply press the SILENT/SET UP button on the top right corner of the base station.


**IF YOU ALREADY HAVE A TELEPHONE PLUGGED INTO THE TELEPHONE JACK AT THE WALL**, unplug the telephone cord from the telephone jack in the wall and plug it into the **GREY** plug on the back of the base station marked **EXTRA PHONE JACK**. Continue with **STEP 5** and **STEP 6**.

**STEP 5**

Plug the telephone cord (provided) into the **YELLOW** jack on the back of the base station marked **PHONE LINE IN FROM WALL OUTLET**. See Fig. 5

**STEP 6**

Plug the other end of the telephone cord into the telephone jack in the wall. See Fig. 6

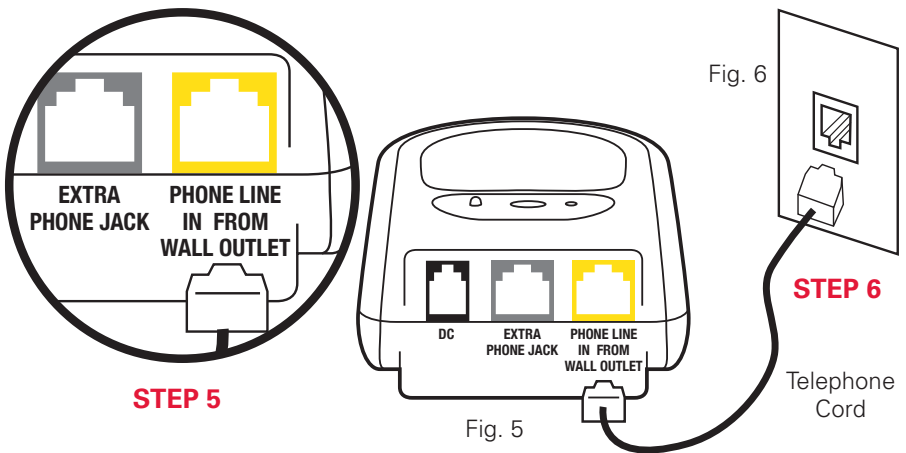If installed properly, the **RED PHONE LINE INDICATOR LED** on the top of the base station should not be illuminated and beeping should cease.

If your telephone system is DSL or VoIP, please see the **MediPendant™ OWNER'S MANUAL** for optional telephone configurations.

**SEE BACK PAGE FOR STEPS 7-10**

## STEP 7

Press and hold the **LARGE ROUND GREY** button for 3 seconds on the front of the pendant. See Fig. 7

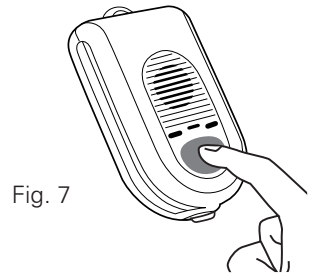
Fig. 7

You will hear the pendant say, **“Your emergency call is now being dialed, please stand by.”** It will say this multiple times. The green light on the pendant will illuminate indicating that it has been activated and will remain illuminated during the entire conversation with the emergency operator. See Fig. 7A

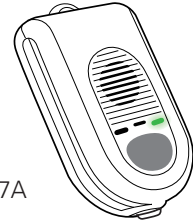
Fig. 7A

## STEP 8

You will then hear a dial tone, followed by a dialing sequence. This is an indication that your pendant is attempting to call the monitoring center.

When dialing sequence is completed, your pendant will again say, **“Your emergency call is now being dialed, please stand by.”**

Soon after, there will be a brief period of silence. This means that your pendant is talking to the monitoring center, and establishing communication.

## STEP 9

Wait a few moments to hear the operator’s voice through the pendant identifying themselves and asking you if everything is OK.

## STEP 10

Identify yourself with the operator and tell him/her that you are simply testing your system. Follow the operator’s instructions to complete the test.

**YOUR MediPendant™ SYSTEM IS NOW ACTIVATED AND READY FOR USE.**

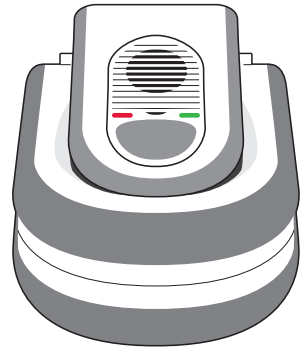
See OWNER’S MANUAL for further instructions.

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